

EN Indications for use: The bubble humidifier is intended to add moisture to breathing gases for administration to patients in homecare, hospital, extended care and hospice environments. Single patient use. Humidity output > 10 mg H₂O/L. Instructions for Use: 1. Add water to the Max level (290 ml). Sterile or distilled water is recommended to reduce possibility of contamination or scaling. If using tap water, boil and cool before using. Replace and tighten the lid and check that it is sealed properly. Change water daily. 2. Attach bubble humidifier to oxygen source. **Important test instructions** refer to graphic 3. To test for leaks: set flow to 3 LPM and block outlet port. Humidifier should sound in 5 seconds or less. If not, check for loose connections. 4. Connect oxygen tubing (not included). If device sounds, check for tubing obstruction. **Cleaning instruction:** Do not sterilize. 5a. Home Use: Every 7 days, wash with mild soap and water, rinse thoroughly. If using boiled tap water in the bubble humidifier, clean every 3 days and air dry. After cleaning the device, visually inspect to ensure device is cleaned. Repeat cleaning steps if necessary. If scaling (mineral deposit) is present, replace humidifier. 5b. Hospital Use: Replace humidifier every 7 days. **WARNING: Leaks may partially or fully reduce oxygen being delivered to the patient.** Caution: Use 3-7 LPM **Contraindications:** None known. For more information: <http://www.salterlabs.com>

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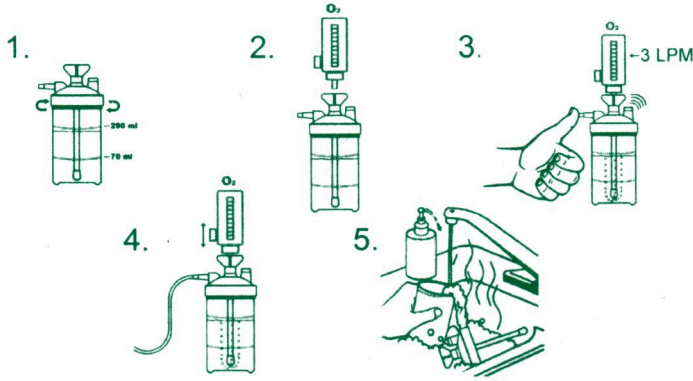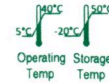(01) 0 0607411 00659 3

LOT (10) 061019A

Manufacturer
SALTER LABS
 30 Spur Drive
 El Paso, TX, 79906 USA
 Made in Mexico

EC REP
 MT Promed Consulting GmbH
 Altenhofstraße 80
 D-66386 St. Ingbert
 Germany

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